

ENROUTE Enflate®

TRANSCAROTID RX BALLOON DILATATION CATHETER



The only transcarotid RX balloon catheter on the US market and the only balloon specifically designed for TCAR


The latest addition to the TCAR® portfolio features:

Short catheter working length

- At 75cm, the shorter catheter working length is optimized for transcarotid access

Balloon lengths optimized for ENROUTE® Transcarotid stents

- Available balloon lengths of 25mm and 35mm are designed for optimal pairing with the ENROUTE Transcarotid stent and minimize unnecessary contact with healthy carotid tissue

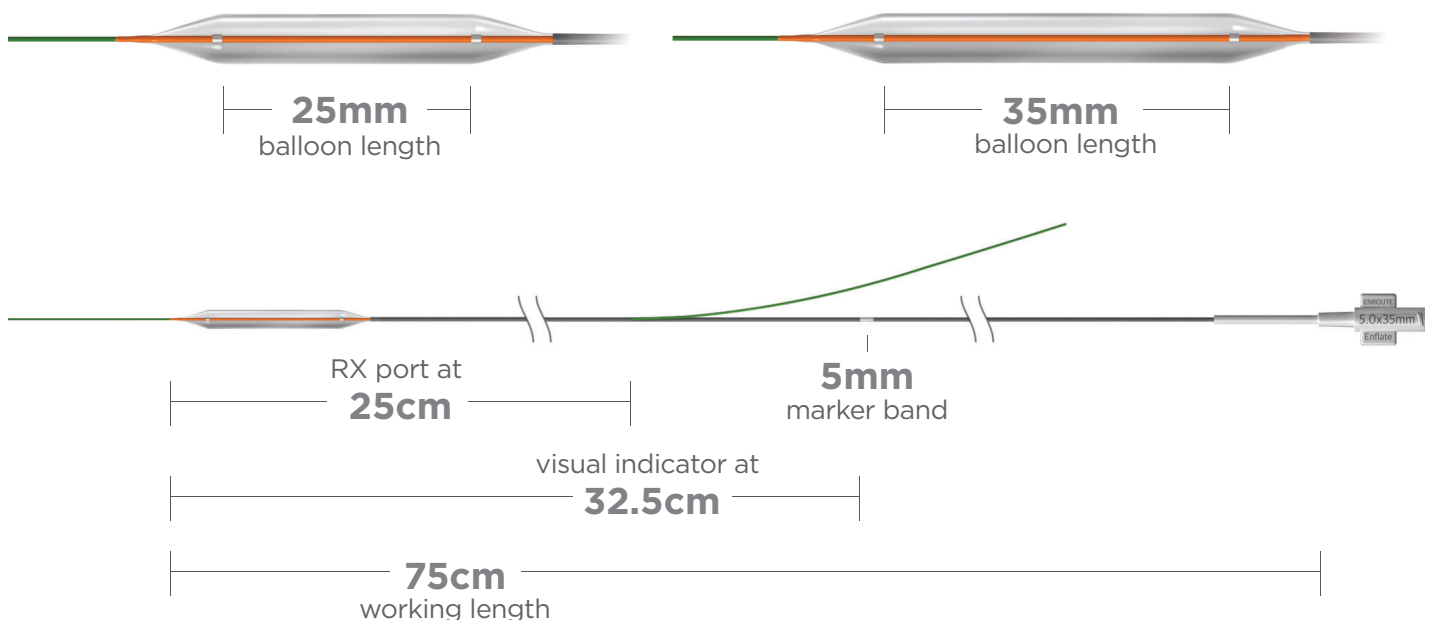

Marker band provides non-fluoro visual indicator that signals when balloon tip is exiting the arterial sheath

- A 5mm-width marker band located at 32.5cm from the tip of the balloon catheter indicates when the balloon is about to exit the arterial sheath during advancement to the lesion
- Designed to indicate when to begin fluoroscopy and minimize fluoroscopy time during angioplasty



Highly visible radiopaque (RO) markers

- RO markers designed to clearly indicate the location of the balloon shoulders for accurate positioning

Additional Features

- Nylon, semi-compliant balloon
- Nominal Pressure: 8 ATM
- Rated Burst Pressure: 14 ATM
- Rapid exchange (RX), 0.014" guidewire compatible
- Low-profile: 4Fr (5Fr for 6.0mm diameter)



Ordering Information: ENROUTE Enflate® Transcarotid RX Dilatation Catheter

Balloon Diameter (mm)	Balloon Length (mm)	Catalog Number
4.0	25	SR-4025-BC
4.5	25	SR-4525-BC
5.0	25	SR-5025-BC
5.5	25	SR-5525-BC
6.0	25	SR-6025-BC
<hr/>		
4.0	35	SR-4035-BC
4.5	35	SR-4535-BC
5.0	35	SR-5035-BC
5.5	35	SR-5535-BC
6.0	35	SR-6035-BC

Indications for Use: The ENROUTE Enflate Transcarotid RX Balloon Dilatation Catheter is intended for percutaneous transluminal angioplasty and post-dilatation of self-expanding stents in the carotid arteries.

Caution: Federal (U.S.) law restricts this device to sale by or on the order of a physician



SILK ROAD MEDICAL, INC. TOLL-FREE: 855.410.TCAR (8227)
 1213 INNSBRUCK DR. MAIN: 408.720.9002
 SUNNYVALE, CA 94089 FAX: 408.720.9013
silkroadmed.com CustomerService@SilkRoadMed.com

